

# Hydrosphere C18



- stable under the use of 100% aqueous eluent
- “hydrophilic” C18 surface for enhanced polar recognition
- no need for ion pair reagents
- based on highly inert, ultrapure, pH neutral silica




Hydrosphere C18	Specification
Particle size / $\mu\text{m}$	3; 5
Pore size / nm	12
Surface area / $\text{m}^2\text{g}^{-1}$	340
Carbon content / %	12
Recommended pH range	2 - 8

## General

The separation of polar compounds in many cases requires highly aqueous mobile phase conditions to achieve sufficient retention on the stationary phase. Conventional reversed phase selectivities do not give reproducible results under these conditions due mainly to the collapse of the C18 chains. Hydrosphere C18 has been developed, on the ultra pure silica support of the ProFamily, as the next generation of speciality phases following the well known YMC-Pack ODS-AQ, which was developed in 1987 and is still a very interesting selectivity option for these purposes.

In order to compare the differences and similarities of both reversed phases, nucleic acid bases are separated under pure aqueous eluent conditions on both selectivities. Interestingly uracil, normally used as a dead volume marker, can be retained on both phases, whereas adenine elutes later on YMC-Pack ODS-AQ. This peak on Hydrosphere C18 is virtually symmetrical with a tailing factor of 1.05, but on YMC-Pack ODS-AQ it shows a slightly tailing peak (with a tailing factor of 1.49).

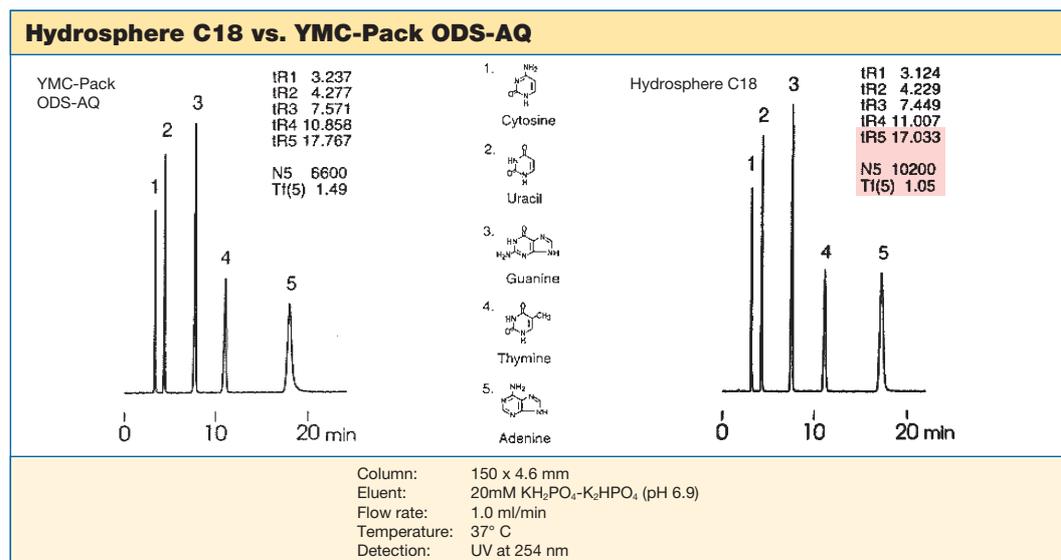

As a consequence, both Hydrosphere C18 and YMC-Pack ODS-AQ are the selectivities of choice within the YMC portfolio to separate polar compounds when using 100% aqueous eluent conditions. YMC-Pack ODS-AQ shows stronger retention for more hydrophobic substances, while Hydrosphere C18 is ideal for basic substances.

Both selectivities are a “must” for approaching the separation of polar compounds.





# Hydrosphere C18

## Properties

A proprietary derivatisation procedure enables Hydrosphere C18 to be penetrated by water without losing its brush-like chain structure. This model of the so-called phase collapse is one explanation for a phenomenon, where the retention times decrease after a very short time under 100% aqueous eluent on a conventional C18 phase.

Here, the retention behaviour of a conventional ODS column and Hydrosphere C18 under 100% aqueous conditions is compared. Its "hydrophilic" C18 surface gives Hydrosphere C18 the capability to show stable retention times even after 24 hours under these chromatographic conditions.

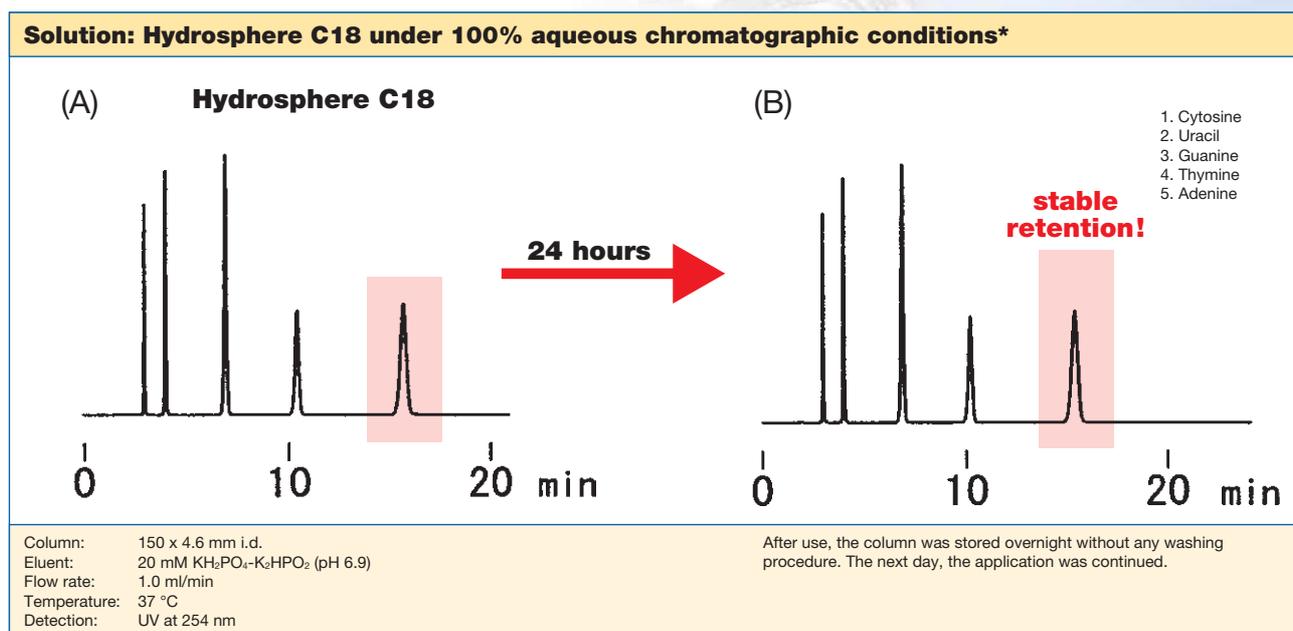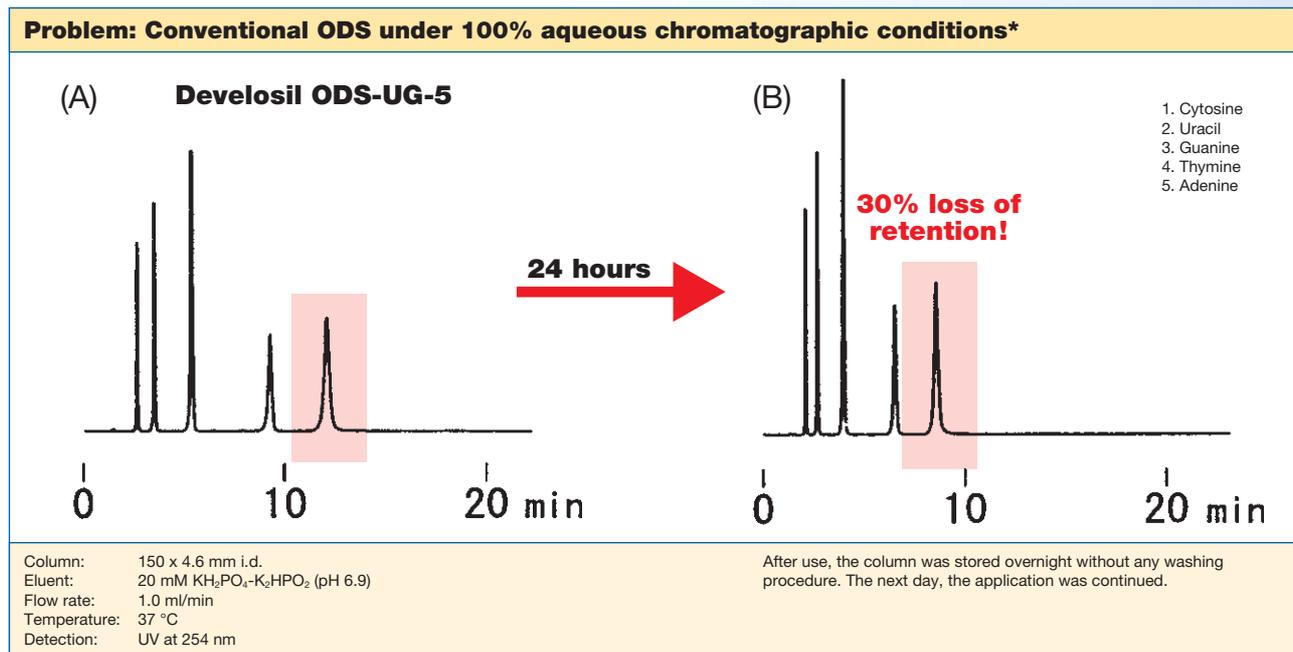




\* By courtesy YMC Co., Ltd.

# Hydrosphere C18

Here, the properties of YMC-Pack ODS-AQ are compared with YMC-Pack *Pro* C18 and Hydrosphere C18 to show similarities and differences of these three selectivities and to give assistance when choosing the right phase for specific chromatographic purposes.



Traditionally, YMC-Pack ODS-AQ was designed for the separation of polar compounds under 100% aqueous conditions and YMC-Pack *Pro* C18 for the separation of acidic and basic substances.

As a consequence, Hydrosphere C18 could be understood as a synergy of the advantages of YMC-Pack ODS-AQ and YMC-Pack *Pro* C18.

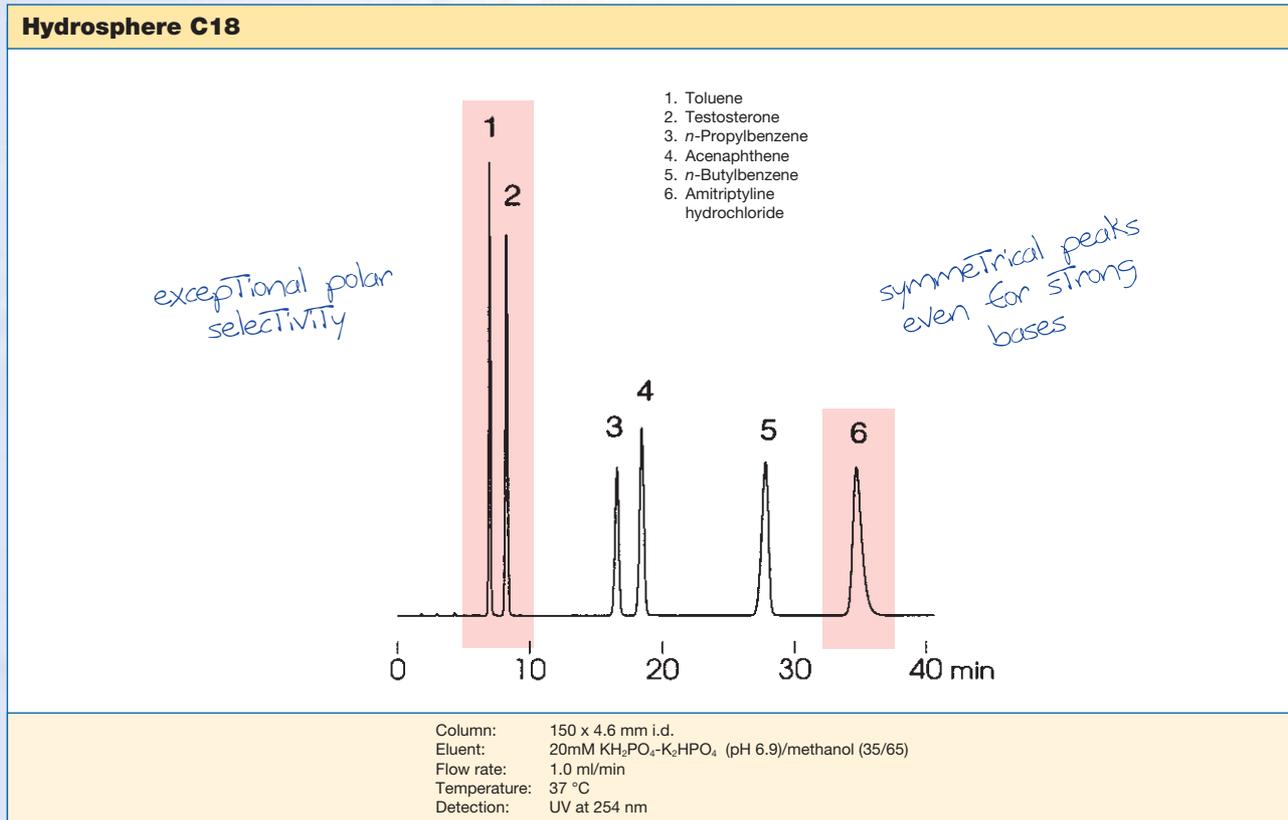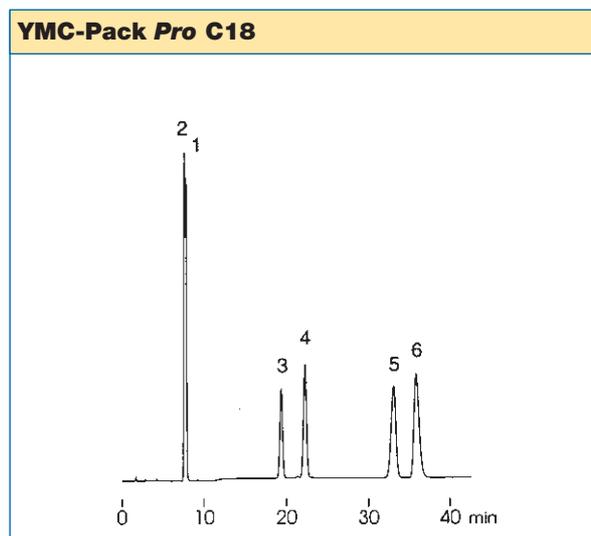
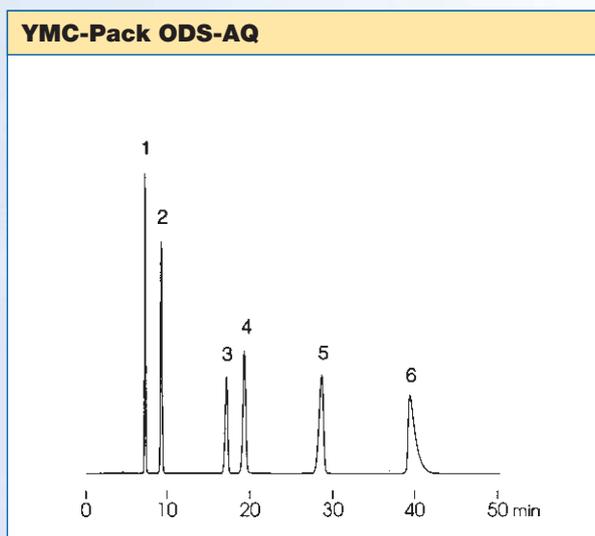


# Hydrosphere C18

Finally, Hydrosphere C18 is, in common with all ProFamily members, applicable to LC-MS methods as shown below. High throughput LC-MS separations with enhanced flow rates are easily achievable with the 3  $\mu$ m version of Hydrosphere C18.

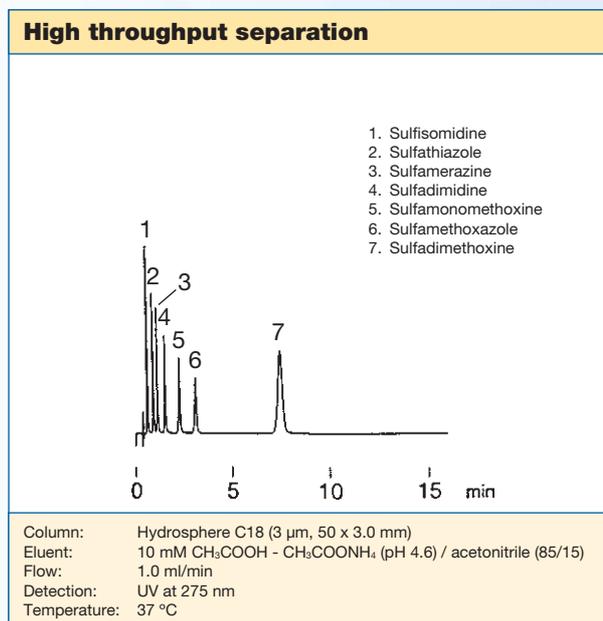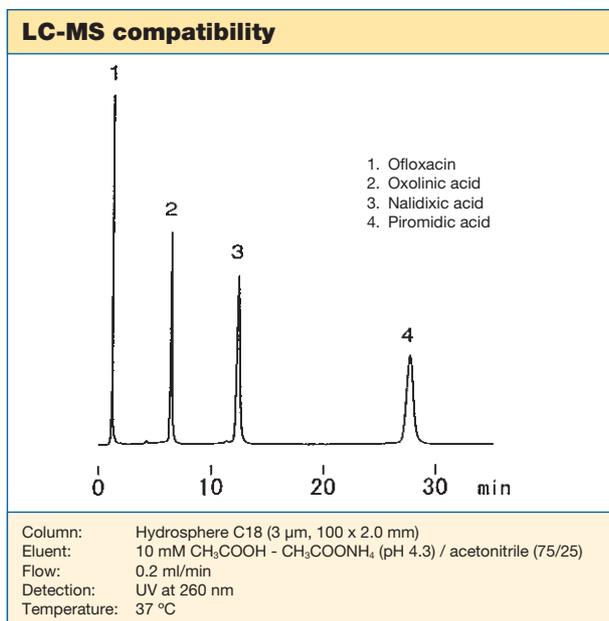

Hydrosphere C18 is the column of choice within the YMC portfolio for a wide range of different substances. As an almost “universal” selectivity it can be used under standard reversed phase conditions as well as pure aqueous eluents as demonstrated by this collection of applications.

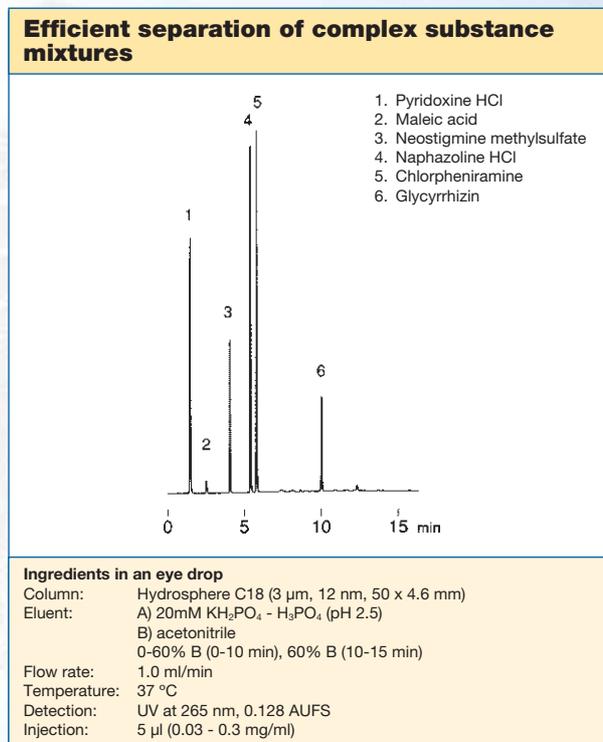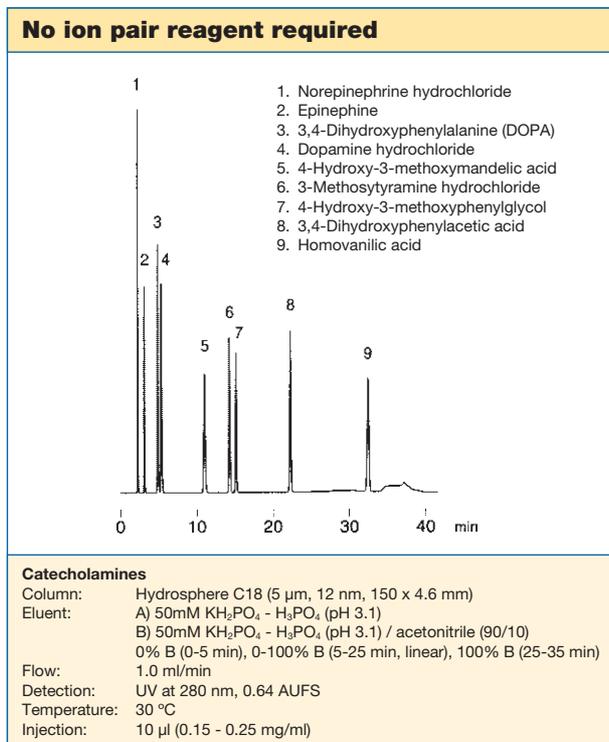

For more applications please refer to our “Application Data Collections” or contact us directly.

## Column care

Hydrosphere C18 is stable towards hydrolysis between pH 2 - 8 in up to 100% aqueous systems and a maximum of 50 °C. Remove acid and buffer salts before storage. Store the column in methanol / water = 70/30.

For detailed information please refer to the “Column Care and Use Instructions” which are shipped with each analytical column.



**High quality & Expert**

# 北京慧德易科技有限责任公司

<http://www.herbs-extract.com>

北京慧德易科技有限责任公司致力于为广大提取物厂家、制药企业、科研机构、大专院校以及化妆品、食品、化工企业提供各种天然提取物和一流的实验室仪器、设备以及试剂、耗材等产品，为客户提供从实验室到工业化生产全过程所需的分离设备和耗材，解决您在分离纯化的生产放大的过程中遇到的问题。

我公司代理的产品有罗门哈斯XAD大孔吸附树脂、三菱化学MCI gel, HP、SP系列大孔吸附树脂、Amersham公司Sephadex LH-20, MERCK色谱产品、TOSOH TSK-gel, Toyopearl凝胶产品, 日本Develosil 和YMC高压液相色谱柱和填料, Sigma-Aldrich标准试剂, 以及Novasep, Applexsion等世界著名色谱设备公司, 我们除了负责以上各公司产品在中国的销售和市场开发以外, 还结合本公司的技术专长, 人才优势, 积极开展分离技术、色谱技术、样品前处理技术的服务工作。

地址: 北京回龙观龙冠大厦719室  
电话: 010-51528296/97/98  
Email: [sales@herbs-extract.com](mailto:sales@herbs-extract.com)

邮编: 102208  
传真: 010-51528299  
<http://www.herbs-extract.com>